




CD-Imager

**The innovative robotic system
for a professional image display
on CD or DVD**

- Fully automatic compact system for creating DICOM patient CDs or DVDs
- Highly compatible with all digital DICOM modalities (multimodality)
- Individual labelling (practice/clinic logo)
- Easy integration of DICOM patient data
- Extremely cost effective due to quick printing times and low ink consumption

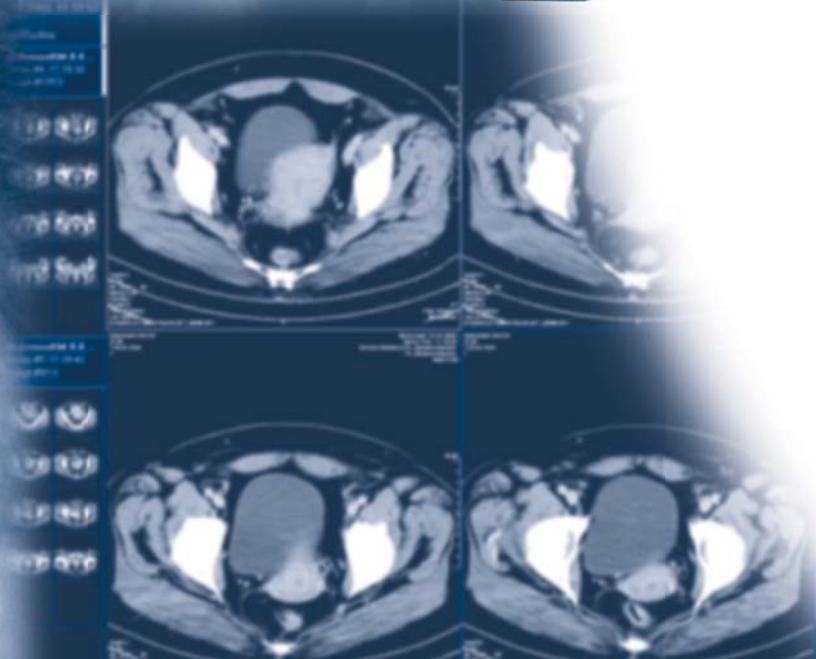

DIGITAL IMAGES ARE OUR WORLD

The medigration CD Imager is a cost effective, fully automatic robotic system for creating CDs or DVDs in DICOM offline media format. It enables you to burn CDs or DVDs of image data and findings for referral physicians or patients in a way that is both

professional yet uncomplicated. The images can be displayed on any PC with the integrated viewer (ImageVision light) in original DICOM quality and can be imported to all PACS systems totally independent from the respective manufacturers.

Extremely easy to use


With our medigration CD Imager it is possible to begin the burning process directly from the modalities or PACS stations within the practice network. As well as patient data, there is the option to transfer the submitted results from the RIS/HIS to HL7. With the easy to manage user interface, you have all running work processes in view all the time and can intervene when necessary. In the sense of corporate identity, our system enables individual labelling by integrating your practice/clinic logo into your print layout. We can incorporate two different logos on request.



Individual labelling

Integrated DICOM viewer

As well as the standard functions (windowing, zoom, inverting, rotation, annotation) our integrated viewer ImageVision light offers the option of displaying automatic cut-lines in images with defined positional relations. In addition, two separate layouts can be set out in the comparison mode. After closing the DICOM viewer (ImageVision light), your desired website opens automatically (homepage of your clinic/practice). On request, your practice/Clinic logo can be integrated.

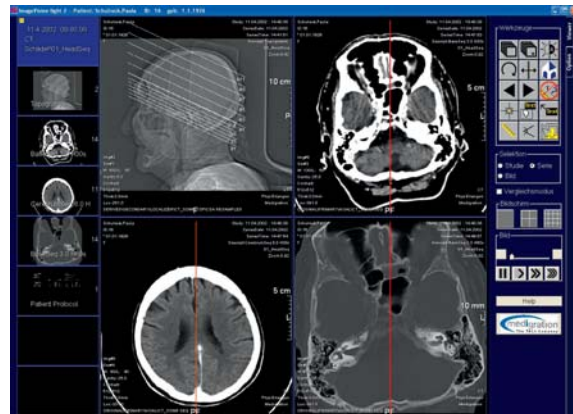

Ergonomic user interface ImageVision

Technical Details

The compact production system has three magazines for 50 media each which allow to finish up to 100 media in one work process and combines ink-jet printing technology, modern burners and transport robotics in a robust, dust-proof casing.

A comfortable drawer eases the taking of single media. The CD Imager is suited for regular use for low to medium production output. Quick variable printing times, low ink consumption and with that, minimal printing costs.

The **money-saving** alternative to laser film!



Gesellschaft für
medizinische Datenverarbeitung
und Integration mbH
Schuhstr. 30
D-91052 Erlangen

Tel: + 49 (0) 91 31 - 6 90 87 - 40
Fax: + 49 (0) 91 31 - 6 90 87 - 50
eMail: info@medigration.de
web: www.medigration.de



For more information please visit
www.medigration.de
If you have any further questions
please feel free to contact us!